

# RC-8DS DISSOLUTION

## RC-8DS DISSOLUTION

RC series of instruments are used for detecting velocity and extent of pharmic dissolution from the tablet, capsule etc in prescriptive menstruum.

### Applicable standards:

National standard( Chinese Codex ) /2015

Corporate standard(Dissolution tester )Q/12XQ0187-2010

### Features:

It has eight vessels and eight poles, with four vessels and four poles in front and other at back.

The head moves up/down automatically.

The basket units and paddle units are made of imported stainless steel (SUS316L).

The magnetism-pumped circular water current can heat the system evenly, and the bath liquid can achieve an equal temperature.


A MPU is used in dissolution tester to automatically control the temperature, rotational speed and time, and the parameters can be preset at any time.

It can be operated conveniently by keystroke. The measurement result is displayed by LED. The preset and real data can be displayed alternately in time.

Automation: auto-test, auto-diagnose, auto-alarm.

It can test by self and alarm automatically. The dissolution parameters of ' Chinese Codex 2010' are stored in and can be taken out at any time.

Time can be preset at many points. The beeper sounds automatically to any point, and then the tester automatically runs to the next one.

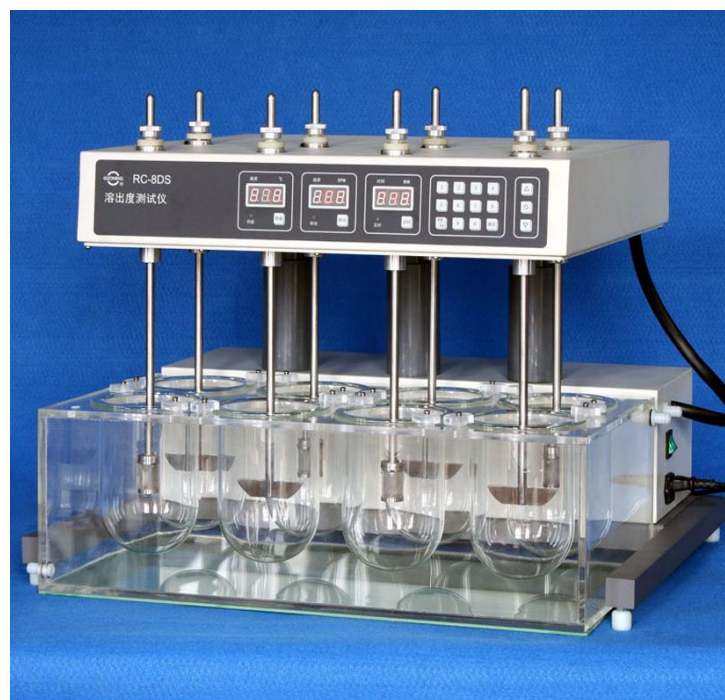

### Specifications:

Speed range: (20~200)rpm

Speed accuracy:  $\pm 2$ rpm

Temperature range ambient to:  $+45.0^{\circ}\text{C}$

Temperature Stability;  $\pm 0.3^{\circ}\text{C}$

Preset timing nine points: (1~999)min

Time accuracy:  $\pm 0.3\%$

Paddle radial run-out:  $\pm 0.5$ mm

Basket radial run-out:  $\pm 1.0$ mm

Power: 220V/50Hz/1200W or 110V/60Hz/1200W

Dimension: (680\*540\*480)mm<sup>3</sup>