


# JY-CZ-BL Rapid SSR vertical electrophoresis tank

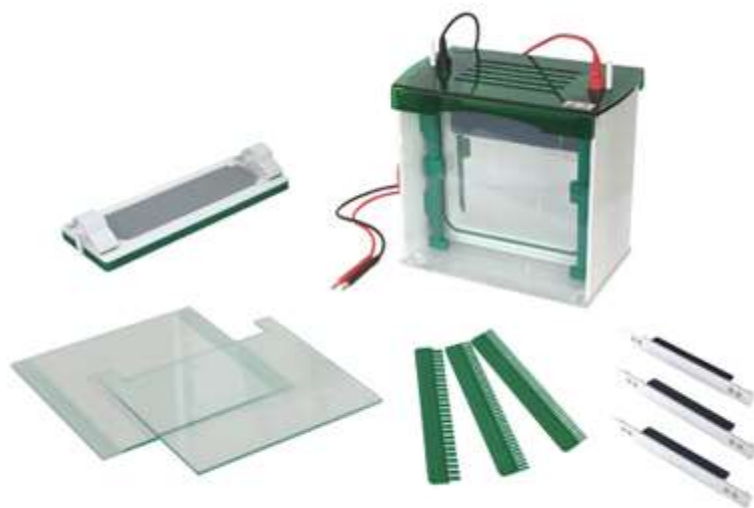

## Technical specifications

Glass plate size (W×L): 216×220 mm

Gel Size (W×L) : 186×205 mm

Gel thickness: 1.0 ( mm )

Gel amount: 1 ~ 2 ( pieces )

Sample volume: ( 1.0mm ) 25 , 40 , 52wells

Buffer Volume: ~ 3500 ( ml )

Dimension ( L×W×H ) : 280×140×250 ( mm )

Net weight: 5.6 ( kg )

## Characteristic features

- The installation of the glass plate in just 15 seconds, and does not need the knob, fast and convenient.
- With the gel making device of normal position, don't should to moving the glass plate from the process of making the gel to electrophoresis.
- With the Precast gel inner core, do not affect the operation of gel making while the electrophoresis.
- Safety cover open button design, convenient opening the cover.
- Abundant buffer not only guarantees cooling effects but also keeps the pH Value stable during the whole experiment process.
- Auto-switch-off when the lid is removed.