


JY-SPAT Horizontal Electrophoresis Cell

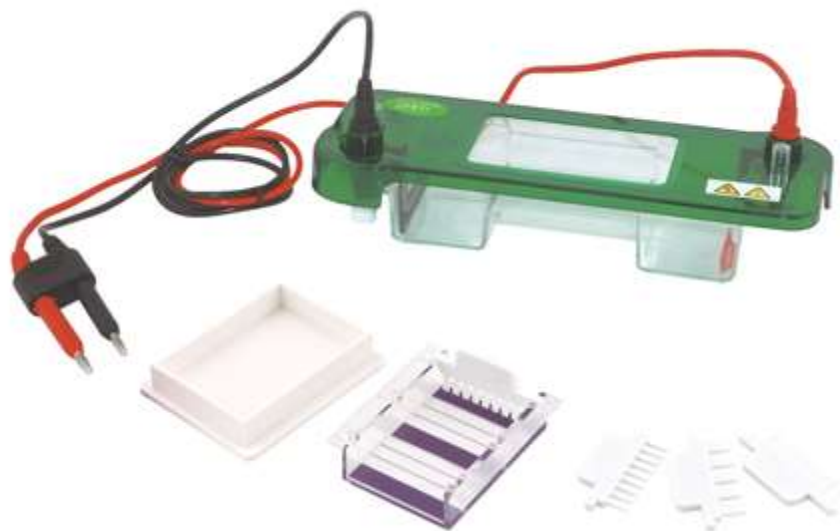

Technical specifications

Gel Size (W×L : 48×75mm

Sample volume: 3, 5, 9.

Buffer Volume: ~150ml.

Dimension (L×W×H: 250×90×50 (mm)

Net weight: 0.4 (kg)

Characteristic feature

- Abundant buffer not only guarantees cooling effects but also keeps the pH Value stable during the whole experiment process.
- Unfolded electrodes make the repair and parts replacement easy and quick.
- With the sample background color design on the tray, convenient to adding the sample.
- The transparent top lid effectively prevents liquid in the cell from volatilizing and avoids electricity leakage.
- Auto-switch-off when the lid is removed.