


JY-SPHT Horizontal Electrophoresis Cell

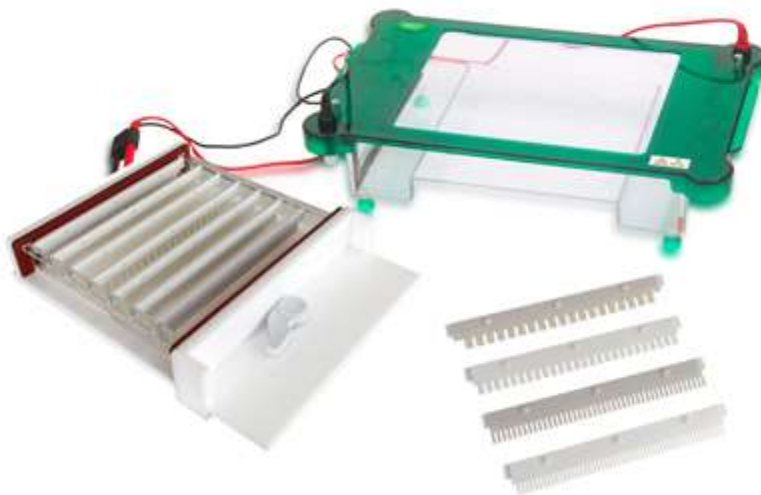

Technical specifications

Gel Size (W×L : 245×250mm.

Sample volume: 19~364

Buffer Volume: ~2000ml.

Dimension (L×W×H) : 487×300×93 (mm)

Net weight: 3.8 (kg)

Characteristic features

- Abundant buffer not only guarantees cooling effects but also keeps the pH Value stable during the whole experiment process.
- Removable electrode enables maintain and replacement easy and quick.
- Adding sample quicker and suitable for the multi-channel pipettes.
- With a special gel casting, preventing the gel leakage, the operation is simple and convenient.
- The transparent top lid effectively prevents liquid in the cell from volatilizing and avoids electricity leakage.
- Auto-switch-off when the lid is removed.